



molift[®]

Molift Mover 180

A light-weight, compact lift, suitable both in home care and institutional environments

Molift Mover 180

A compact, lightweight mobile lift with great flexibility and large lifting range. Suitable in all care environments, from homes to hospitals.

Large hoisting range

Molift Mover 180 has a lightweight aluminum construction that makes the lift very easy to maneuver. The low weight in relation to lifting capacity is best in class. A large lifting range allows lifting from high surfaces as well as from the floor. The integrated Molift Service Tool ensures correct service and maintenance intervals.

Ergonomic handling

For improved caregiver ergonomics, the Molift Mover 180 is equipped with Molift's unique intuitive push handle. The handle provides multiple and flexible grip possibilities for improved caregiver ergonomics. The easily accessible hand control is very user friendly.



Watch our videos for inspiration and instructions

| Product | Item no. | Dimensions | Technical |
|--|--|---|---|
| Molift Mover 180 - with 4-point sling bar and 4 inch castors | M18003 | Hoisting height with 4-point sling bar Outside measurement with closed legs Outside measurement with open legs | NiMH 14,4 V – 2,6 Ah 3 hours 1.81 inch/sec 50 lifts/165 lbs /19.7 inch |
| Material Aluminium, steel, plastic | | Technical Battery Battery charge time Lifting speed 165 lbs Rated performance | |
| Weight Safe working load (SWL) Weight of unit (excl. sling bar and battery) | 400 lbs 70.5 lbs | Cleaning | |
| Dimensions Height Length Wheelbase width Leg height Castor size Hoisting interval Hoisting height with 2-point sling bar | 77.6 inch 49.6 inch 26.4 inch 5.3 inch 4 inch 53.8 inch 14.6 - 68.1 inch | | |
| | | Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution | |

Features



Ergonomic push bar

The unique push bar provides the best possible natural grip, creating good ergonomics and many grip possibilities.



Sling bar options

The lift can be used with both a 2-point and a 4-point sling bar. The sling bars are available in many variants and can easily be replaced without tools.



Molift Service Tool

Built-in software which calculates the number of lifts and gives a warning signal when service is required. It can also be used to retrieve data about lifting usage and lifting patterns.



Easy to manoeuvre

Molift Mover 180 has a total weight of only 70.5 lbs which, together with the design, makes it easy to manoeuvre.




Molift weight scale

The Molift weight scale is approved Class III and complies with all technical standards. It can also calculate BMI (Body Mass Index) with a simple touch.



Sling bars

The Molift Mover 180 can be used with 2-point and 4-point sling bars, available in a variety of sizes. The quick-release coupling makes it easy for the carer to easily change the sling bar.

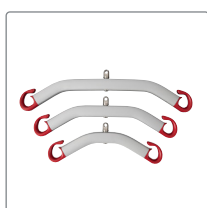

2-point sling bars

| Size | Item no. | W x H (inch) |
|------|----------|--------------|
| S | 1830043 | 13.8 x 5.1 |
| M | 1830042 | 17.7 x 5.1 |
| L | 1830041 | 21.7 x 5.1 |



4-point sling bars




| Size | Item no. | L x W x H (inch) |
|------|----------|-------------------|
| S | 1830110 | 13.8 x 11.8 x 8.3 |
| M | 1830111 | 17.7 x 11.8 x 9.5 |
| L | 1830112 | 21.7 x 13.7 x 9.5 |
| XL | 1830113 | 25.6 x 14.6 x 9.5 |



Etac is a world-leading developer of innovative assistive devices and patient handling equipment. Since 1973, we have been committed to improving quality of life for the individual, family members and caregivers.

For the intended purpose of the products and latest information, please visit www.etac.com

Etac

-  +1 (833) 818-0024
-  sales.patienthandling@etac.com
-  www.etac.us.com


Creating Possibilities